

Avon 120™

Previously Norfolk Bari™

The Avon 120™ Bariatric Mattress has been specifically designed to suit the differing weight and body shapes of bariatric patients. The mattress profile design has been developed to support the heavier patient group while optimising pressure re-distribution.


The multi-stretch, vapour permeable, high frequency welded cover ensures compliance with strict infection control policies, whilst safeguarding the skin integrity of the patient through a reduction in shear and friction forces.

The Avon 120™ has been designed for those patients considered to be at 'High Risk' of pressure ulcer development.

**If you have any questions we are here to help.
Call us on 0800 193 6030.**

PRODUCT FEATURES

- Designed to suit the differing weight and body shapes of bariatric patients
- Developed to support heavier patients while optimising pressure re-distribution
- Multi-stretch, vapour permeable, high frequency welded cover

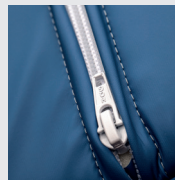

Features and Benefits

- Multi-stretch, vapour permeable, high frequency welded cover with non-slip base (sewn cover also available). Aids infection control, reduces shear & friction and improves patient safety.
- Covers can be easily cleaned and laundered up to 80C.
- Cut foam surface ensures total comfort is achieved, whilst delivering excellent pressure re-distribution.
- Fire retardant CMS foam fully complies with British Standards.
- UK Manufactured product. Provides total confidence to user via exceptional quality and workmanship.
- Clinically proven. Giving peace of mind to nursing staff and patients alike.
- 50-V Insert.

TECHNICAL INFORMATION

| | |
|---------------------------|--|
| Product Code | AV120 |
| Risk Category | High Risk |
| Fatigue Class | Very Severe |
| Maximum Patient Weight | 50 stone 317kg |
| Mattress Dimensions | 198cm x 120cm x 15cm (Bespoke/made to measure service available) |
| Mattress Weight | 18kg |
| Foam Density | Insert: 38 – 40kg / m ³ |
| Nominal Hardness | Castellated Insert: 185 – 210N |
| Mattress Cover Coating | 205 – 235g (choice) |
| Fire Retardancy Standards | Exceeds minimum requirements of Crib 5 Composite testing (certification available upon request) BS7177 |
| Warranty | 2 years – extended warranty available |