


X'PRESS®



Bed with MADELIA boards.
Light Beech finish and epoxy side rails and lifting pole

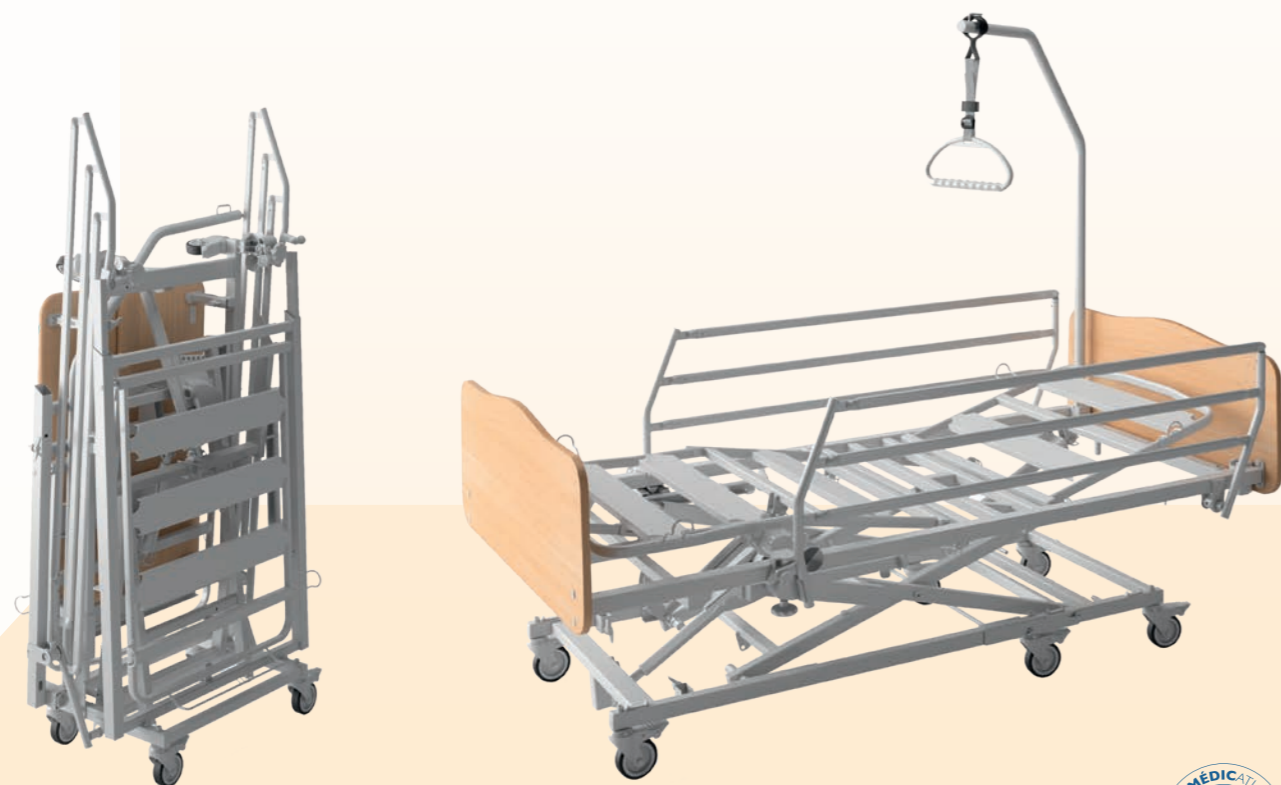

Transformation of the self-supporting trolley into a frame.



The foot part of the bed base are locked in horizontal position.



The two parts of the bed base is joined securely together with through rondo screws.



Boards are fastened using the **EASY MOVE** system: easy, fast and secure.



transto X'press kit flat loading.

ERGONOMICS

- ▶ Thanks to the bed's mechanical ergonomics the user is able to adopt the recommended postures to help prevent MSDs (Musculoskeletal Disorders) when the bed is being delivered, assembled and dismantled.
- ▶ Because the bed does not take up much space and is stored on its own castors, it is easy to move in any situation (elevators, small vehicles, etc.).
- ▶ This self-supported bed does not need any trolley.

EASE OF USE

- ▶ Assembly and dismantling operations are very implicit. No tools are required.
- ▶ Blue parts play a functional role. They provide guidance for the bed's assembly.

SAFETY

- ▶ The new Easy Move board fastening system avoids loss of any rondo screws and enables the bed boards to be fastened securely.
- ▶ The cable routing reduces pulling and risks of deterioration.

DESCRIPTION	DIMENSIONS	ADJUSTABLE HEIGHT	BACK REST	LEG REST	LIFTING CAPACITY	CASTORS	STANDARDS	MOTORIZATION
X'PRESS 2 functions	Assembled bed: 90 x 200 cm divisible bed base	34 to 84 cm - Ground clearance: 13 cm	0° to 75°	12° manual cranks	SWL ⁽¹⁾ 170 kg mattress 20 kg patient 135 kg accessories 15 kg	Ø 100 mm	IEC 60601-2-52 compatible with environments 3, 4 and 5	LINAK IP66
X'PRESS 3 functions without knee break	Bed in storage position: 50 cm wide, 90 cm long and 160 cm high			13.5° electric				
X'PRESS 3 functions with knee break				thigh ection 26° and 16° feet adjustment by manual crank				

⁽¹⁾ Safe working load

- 2 accessory holders (lengthwise)
- 2 IV stand holders (at the head of the bed)
- 2 lifting pole (at the head of the bed)

Handset, electric and mechanical parts: 5-year warranty*

Boards and accessories: 1-year warranty*

*Under normal conditions of use and excluding wearing parts - please see instructions of use

Photos not contractually binding

FP000001/01.2018

